



Cryo Multimaster 6 in 1 skincare device

Electroporation



Cooling Therapy



Heating Massage



Sonophoresis



Iontophoresis



Skin Scrubber



MULTI MASTER

Cryo Multimaster is that aesthetic medical device capable of delivering all core functions required for various skincare program . It has 6 set of functions.

1. **Cryo-Electroporation Care:** It is called electric puncture method. This function enables drug penetration to the desired site by using galvanic current without using a needle. It applies electrical pulses to the skin surface to disturb the skin's dual phospholipid membranes, temporarily creating holes to penetrate the drugs needed for the skin regeneration.
2. **Cryo-Cooling therapy Management:** It is applied immediately after the procedure that causes thermal reaction by freezing electrophoresis treatment. It reduces the erythema and the swollen sensation thereby calming the swollen skin. Various regen ampoules can be applied for faster skin regeneration.
3. **Cryo- Heating massage Management:** It delivers muscle relaxation, improves skin elasticity, pore relaxation, stimulate blood circulation and promote metabolism.
4. **Sonophoresis:** Uses ultrasonic vibrations to deliver nutrients and moisturizing ingredients to deep layers of skin.
5. **Iontophoresis:** Penetrate ionized vitamin C, collagen, ceramide solution into the skin cells using galvanic current of electro-ion migration device. This method treats various pigmentations that has been difficult to threat through other techniques.
6. **Skin scrubber:** The vibration of ultrasonic wave extracts waste and fine dust deposited inside pores and skin.

MULTI MASTER


1. Cryo Electroporation

AKA electric puncture method. It penetrates drug to the desired site by using galvanic current without needle. It applies electrical pulses to the skin surface, break the skin's dual phospholipid membranes, temporarily creating holes to penetrate drug needed for skin regeneration.

Principle	Boosts medicine absorption by making instant holes on the skin with Electric pulses.
Effect	Increases medicine permeation, soothe skin by the cooling, promotes Metabolism by the heating.
Benefits	Boosts drugs delivery by 70%.
Method	Rub with cyro after applying medicine on the skin.



Cryo

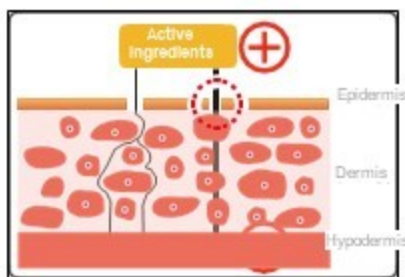

Electro poration


MULTI MASTER

2. Cryo Cooling Therapy Management

It is applied immediately after the procedure that causes thermal reaction by freezing electrophoresis treatment. It reduces the erythema and the swollen sensation thereby calming the swollen skin. Various regen ampoules can be applied for faster skin regeneration



Cryo

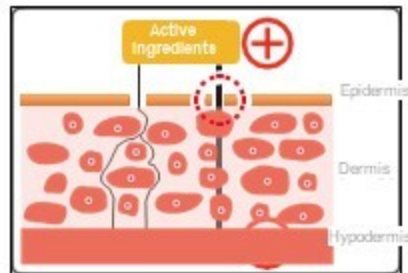

Electro poration



Cooling Therapy


MULTI MASTER

3. Cryo- Heating massage Management

It delivers muscle relaxation, improves skin elasticity, pore relaxation, stimulate blood circulation and promote metabolism.



Cryo

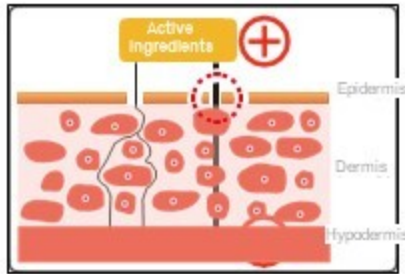

Electro poration



Cooling Therapy



Heating Massage

MULTI MASTER

4. Sonophoresis - Ultrasonic Wave Care

Discharges dust and other skin wastes by ultrasound wave vibrations and helps make clear skin.

Principle	Medicine permeation by the ultrasound wave vibration, promotes cell restoration.
Effect	Skincare by heating, vibration, cleansing, permeation.
Benefits	Improves Wrinkles, increases skin elasticity, Fat care
Caution	Be careful use on the scars, tumors and around bones.
Method	Use after applying the ultrasound transmission gel on the spot.



1. 24,000Hz-28,000hz the face care, dead skin cell removal
2. 1mhz : face and body care, more on the body care Fat body care – thigh, hip, calf and belly effective on the highly dense fat areas
3. 3mhz : good for deep cleansing and it gently massages hypodermis with the scrubber vibrating 3,000,000 times per second to the 1.5cm deep inside of the skin. It stimulates the collagen and elastin layer and helps blood and lymph circulation. This way, it helps skin revitalization and anti aging care from the cell restoration.


MULTI MASTER

5. Iontophoresis

Dark Spot and melanin pigmentation care by penetrating the ionized vitamin c on the spot using the galvanic electro action of the iontophoresis device.

Principle	Penetrate the ionized medicine into the skin by the galvanic electro action.
Effect	Whitening, healing skin
Benefits	No pain, easy to use, good for the spot care
Method	Proceed with applying ionized material on the spot.

Iontophoresis

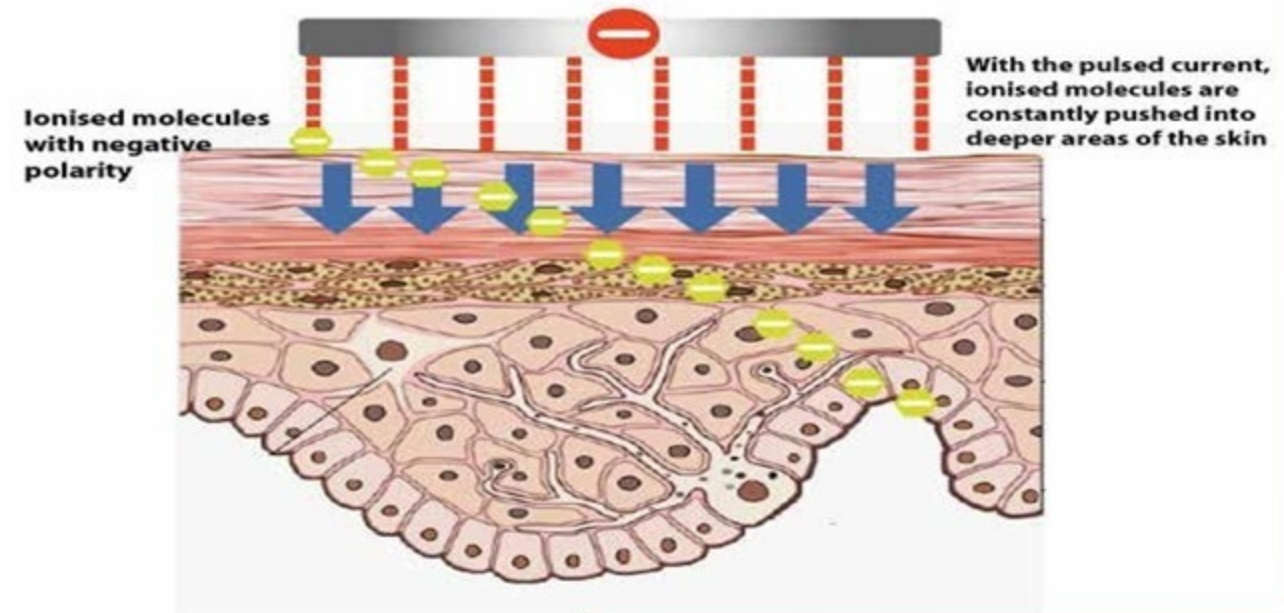

MULTI MASTER

6. Scrubber- Ultrasonic Wave Care

Discharges dust and other skin wastes by ultrasound wave vibrations and helps make clear skin.

Principle	Discharges the remaining dusts and skin wastes by ultrasound wave vibrations.
Effect	Helps keep the skin clear.
Benefits	No pain, easy to use, good for the spot care
Method	Use after applying the ultrasound transmission gel on the spot.



	Vibration	Scrubber edge usage
M1 : Dead skin cell removal		
M2 : Nutrition		
M3 : Wrinkle care		